



IRIDEX

Elegantly Simple Solutions™

# Why 577 nm Yellow?

Clinical Benefits of 577 nm Yellow Laser  
in the Treatment of Ocular Disorders

IQ 577™ Laser Benefits

Overview of Ocular Photocoagulation Lasers  
with Yellow Wavelengths

Understanding Laser-Tissue Interaction, Absorption,  
and Conversion of Laser Energy into Heat

MicroPulse™ Technology for  
Tissue-Sparing Photocoagulation

Laser Treatment Strategies




# IQ 577™ Laser Benefits

1. High transmission through dense ocular media<sup>1,2</sup>
  - Longer wavelength with lower light scattering results in less power needed for the intended retinal irradiance
2. Consistent laser lesions for fast procedure time (See Figure 3)
  - Consistent tissue uptake and reduced thermal spread
  - Less frequent need to readjust laser parameters
3. Enhanced visibility for reduced intraretinal damage<sup>2</sup>
  - Enables early observation of very light tissue reactions at the level of the retinal pigment epithelium (RPE)
4. Low power required for increased patient comfort<sup>9</sup>
  - Lower transmission to deeper tissues<sup>2,4</sup>
5. Allows treatment closer to the macula
  - Negligible absorption by xanthophyll<sup>2</sup>
6. Most efficient focal treatment of vascular structures (See Figure 1)
7. Tissue-sparing capability through MicroPulse™ technology



Figure 1. Laser Wavelength & Effective Light Absorption

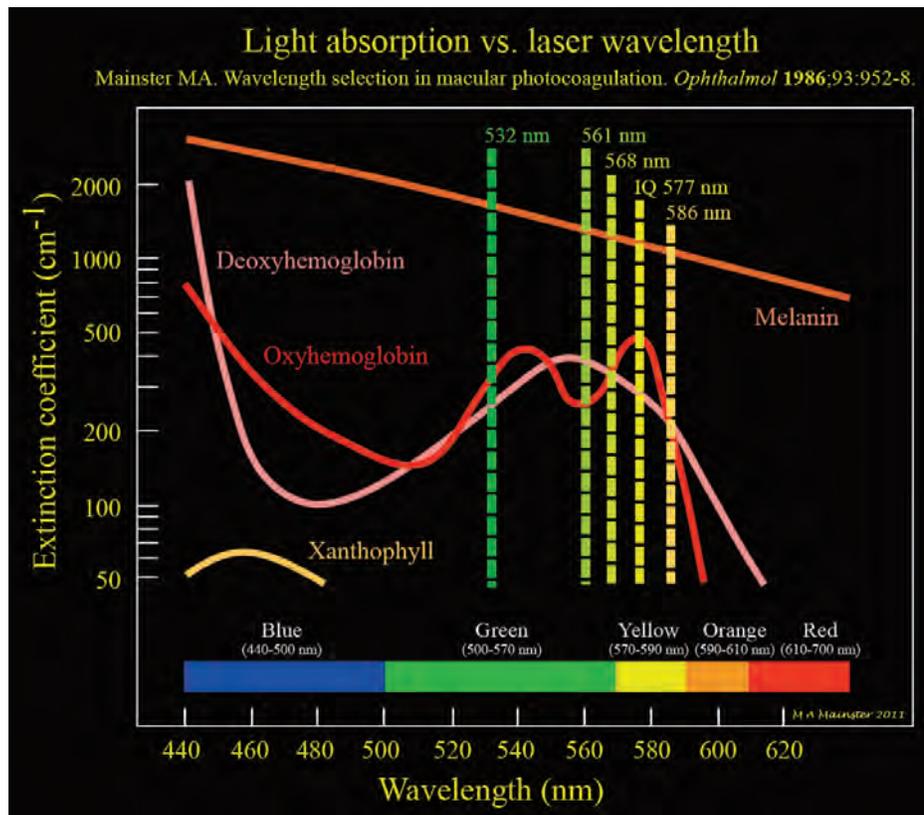

Laser	$\lambda$ (nm)	Extinction Coefficient cm <sup>-1</sup>			Ratio HbO/HbR	Ratio HbO/Melanin
		HbO	HbR	Melanin		
Argon Green	514	150	160	1850	0.94	0.08
FD Nd:YAG Green	532	320	250	1600	1.28	0.20
DPSS "Yellow"	561	250	375	1300	0.67	0.19
Krypton "Yellow"	568	330	330	1200	1.00	0.28
<b>IRIDEX Yellow</b>	<b>577</b>	<b>460*</b>	<b>275</b>	<b>1130</b>	<b>1.67</b>	<b>0.41<sup>^</sup></b>
FD Nd:GdVO <sub>4</sub> Yellow	586	210	210	1040	1.00	0.20

\* 577 nm is at the absorption peak of HbO, which is important for direct treatment of vascular structures.<sup>2</sup>

<sup>^</sup> 577 nm has the highest ratio of HbO to melanin extinction to minimize damage to underlying pigmented tissue.<sup>2</sup>

## Overview of Ocular Photocoagulation Lasers with Yellow Wavelengths

Argon-dye laser photocoagulators were introduced in the mid-1980s and quickly became popular with retina specialists. Using argon-dye photocoagulators, clinicians could select green (514 nm), and a tunable range of wavelengths providing greenish-yellow to orange to red (560 nm to 630 nm); among them, the 577 nm (yellow) - at the peak of the oxyhemoglobin absorption curve - quickly became the favorite wavelength for a variety of reasons including less scatter and increased efficiency. Eventually, argon-dye photocoagulators became less popular because they were costly, complex, and difficult to maintain.

The next yellow laser wavelengths became available in a 3-color krypton laser that delivers 568 nm, followed by solid-state lasers that deliver 561 nm, 568 nm, 577 nm, and 586 nm – all marketed as “yellow” light; however, not all yellow wavelengths are alike: 561 nm and 568 nm are in the *green* wavelength spectrum (500 nm to 570 nm). Although both 577 nm and 586 nm are in the yellow wavelength spectrum (570 nm to 590 nm), 577 nm has higher absorption coefficients in oxyhemoglobin (HbO), deoxyhemoglobin (HbR), and melanin. (See Page 3)

The importance of laser wavelength absorption characteristics is discussed below.

## Understanding Laser-Tissue Interaction, Absorption, and Conversion of Laser Energy into Heat

1. There are three principal chorioretinal light-absorbing chromophores:

### Melanin

Melanin is the most effective light-absorbing chromophore. It's located in the RPE and choroid, where light energy converts into heat. Light absorption in melanin decreases with increasing wavelengths. Melanin concentration varies among patients and fundus locations, producing variability in light absorption. (See Figure 2)

### Hemoglobin


HbO and HbR are important absorbers after melanin. Their absorption spectrum (See Figure 1) is characterized by distinctive peaks: 542 nm green and 577 nm yellow in HbO; and 555 nm in HbR. High choriocapillaris hemoglobin absorption provides more uniform laser effects in patients with light or irregular fundus pigmentation. (See Figure 3)

### Xanthophyll

Xanthophyll is located in the inner and outer plexiform layer of the macula where thermal damage is undesirable. 577 nm is minimally absorbed by xanthophyll, so there is negligible light absorption in the inner retina or its resultant temperature elevation.

2. Laser energy is absorbed in the RPE from which heat can spread to the overlying neurosensory retina. When the retina becomes thermally damaged it loses its transparency and scatters white slitlamp light back at the observer which appears as a “blanching” endpoint. Higher temperature results in greater loss of retinal transparency with increased scattering and whiter endpoint. (See Figure 4)

Figure 2.  
Variability of Melanin Distribution

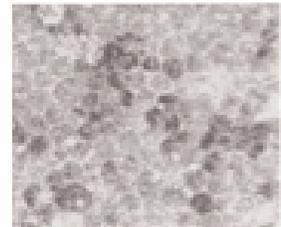


*Human RPE shows irregular distribution of melanin.*  
Geeraets WJ, et al.  
The relative absorption of thermal energy in retina and choroid.  
Invest Ophthalmol  
1962;1:340-7

**Figure 3. Melanin and HbO Absorbing Chromophores Increase Efficiency of 577 nm**

577 nm has the highest combined absorption in the melanin-oxyhemoglobin layers of the RPE/choriocapillaris complex. (See Figure 1)


*Melanin is unevenly distributed in the RPE and choroid. Based on the laser wavelength, some light will absorb, and some light will pass through.*

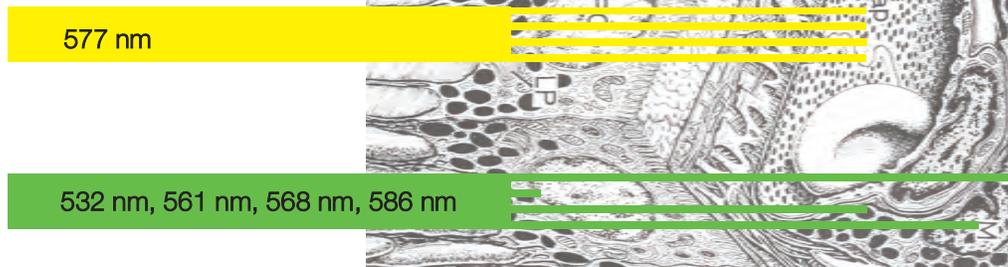

*Hemoglobin in the choriocapillaris is more uniformly distributed for a more consistent uptake of laser light. 577 nm has the highest absorption coefficient in HbO. (See Figure 1)*

*The lower absorption and increased transmission of 577 nm through the non-uniform melanin granules of the RPE is more than compensated by the higher absorption of 577 nm in the underlying more uniformly distributed hemoglobin-rich choriocapillaris.*

**Figure 4. Conversion of Light Energy into Heat**

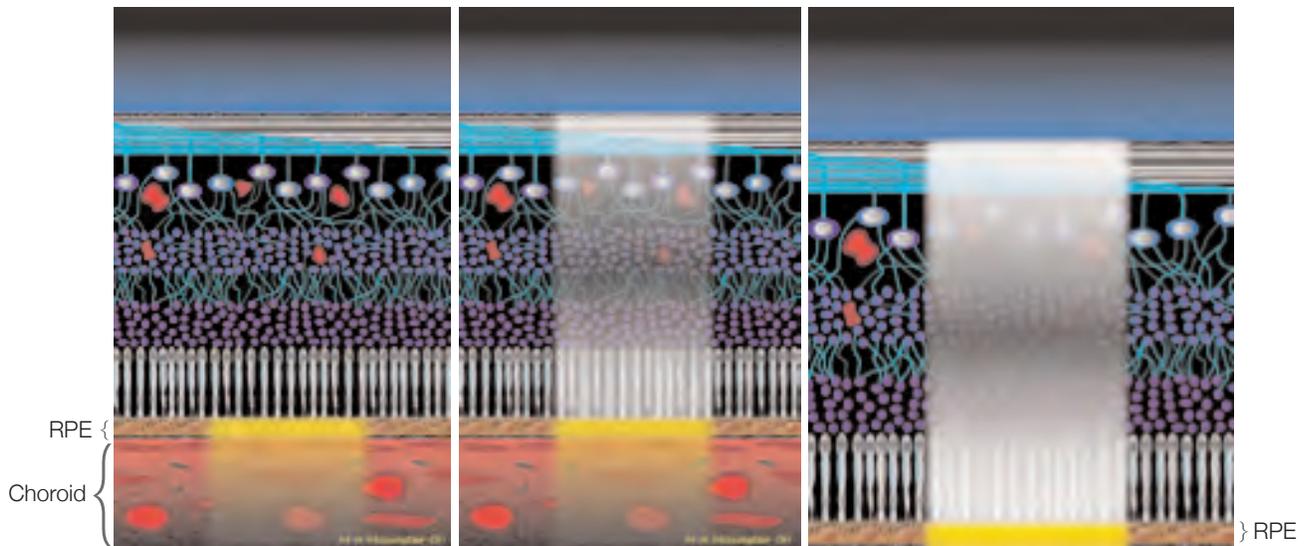

*Heat conduction spreads temperature rise from laser light-absorbing pigmented tissues (melanin in the RPE and choroid, and hemoglobin in the choriocapillaris in the choroid) to overlying neurosensory or collateral retina.*

*The overlying retina damaged by heat conduction loses its transparency and scatters white slit-lamp light back at observers.*


*More damage means less transparency and a whiter lesion.*

*Illustrations compliments of Martin A. Mainster, PhD., MD, FRCOphth*

# MicroPulse™ Technology for Tissue-Sparing Photocoagulation

MicroPulse is a tissue-sparing laser technology and dosing protocol that can limit thermal elevation to temperatures below the threshold of retinal tissue damage to induce beneficial intracellular biological effects without any visible laser-induced damage during and at any time post treatment. (See Figure 5) Early MicroPulse protocols used 810 nm, and have been refined over the years to further improve treatment outcomes while continuing to do no harm. More recently, studies using 577 nm MicroPulse protocols have been presented.

Figure 5. MicroPulse Technology

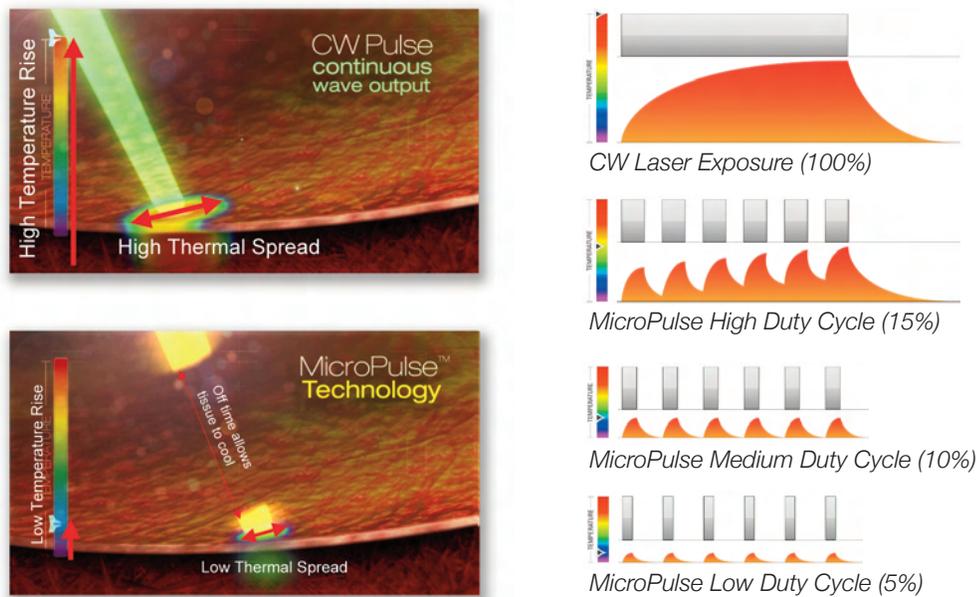

*In conventional, continuous-wave (CW) photocoagulation, a rapid temperature rise in the target tissue creates blanching and a high thermal spread. MicroPulse technology finely controls thermal elevation by “chopping” a CW beam into a train of repetitive short pulses allowing tissue to cool between pulses and reduce thermal buildup.*

The clinical efficacy of MicroPulse protocols has shown favorable therapeutic responses with minimized collateral effects in the treatment of diabetic macular edema (DME),<sup>5-19</sup> proliferative diabetic retinopathy,<sup>6, 20, 21</sup> macular edema secondary to branch retinal vein occlusion (BRVO),<sup>6, 22-24</sup> and central serous chorioretinopathy (CSC).<sup>25-30</sup>

In a prospective, randomized, controlled clinical trial for DME, MicroPulse photocoagulation has demonstrated to be as effective as conventional (Early Treatment of Diabetic Retinopathy Study) photocoagulation in **stabilizing visual acuity and in reducing macular edema, with the benefits of no tissue damage detectable at any time point postoperatively and of significant improvement in retinal sensitivity.**<sup>15</sup>

In 2010, the first study was reported on the use of the IQ 577™ in its MicroPulse mode to treat CSC. Results showed all patients had complete resolution of their symptoms, there was functional visual improvement, and there were no signs of laser marks on the treated areas by clinical examination or fluorescein angiography.<sup>31</sup> Additional studies using 577 nm tissue-sparing photocoagulation have shown clinical effectiveness for the treatment of CSC,<sup>32</sup> DME,<sup>33, 34</sup> and BRVO.<sup>35</sup> (See Figure 6)

# Laser Treatment Strategies

## 577 nm Continuous-Wave Mode

The 577 nm wavelength photocoagulator has been shown to produce visible endpoints similar in both appearance and clinical efficacy to those achieved using green laser photocoagulators, while requiring only 60% to 70% of the green laser power.<sup>1,3</sup> Therefore, a useful way to familiarize yourself with the ophthalmic response to 577 nm wavelength strategy might be to initiate each laser treatment with 50% to 60% of the power you would normally use with your green laser; keep all other techniques and parameters the same. Titrate power to achieve your desired endpoint. Note that focal treatment of microaneurysms may be accomplished using very low powers (90 to 130 mW) compared to other visible wavelengths.<sup>36</sup> Anecdotal experience has shown powers as low as 50 mW.

*Note: The IQ 577™ laser has a user preference selection to automatically turn off the aiming beam during treatment laser emission. This allows you to more easily observe the earliest and most subtle visible tissue reactions at the level of the RPE without the distraction and contrast reduction due to scattered aiming beam light.*

## MicroPulse™ Mode


MicroPulse is typically used to administer subvisible threshold laser treatments to macular and perimacular targets. When used here, the terms “subvisible,” “subvisible threshold,” or “subthreshold” denote that the desired endpoint is one in which treated tissue offers no ophthalmoscopically observable laser effects. Nevertheless, 577 nm and 810 nm studies have confirmed that subvisible laser treatment strategies can be clinically effective while inducing no tissue changes discernable by slitlamp observation, fluorescein angiography (FA), or fundus autofluorescence (FAF) at any time post-operatively.<sup>15, 35</sup> Subvisible MicroPulse laser treatments are consistently effective without causing such changes because the total laser energy is only a percentage (often chosen by clinicians to be 20-70%) of that needed to produce a visible endpoint.

Energy (J) is equal to [Laser Power (W)] x [Exposure Duration (s)] x [Duty Factor (%/100)]. Duty Factor is often 5% to 15% when using MicroPulse mode, and is 100% when using CW mode. Clinicians have reported various strategies to adjust these parameters relative to suprathreshold burns in order to achieve clinically effective subvisible endpoints.<sup>13, 15, 19, 35</sup>

Additional parameters to consider in any laser treatment protocol, and particularly during MicroPulse, is spacing between laser treatment spots, and the total number of treatment spots administered. Due to the limited thermal spread of MicroPulse exposures, subvisible treatments often call for the administration of a greater number of treatment spots with denser spacing than used for threshold laser grid treatments.<sup>17</sup>

For more information on MicroPulse, register at [www.iredex.com/micropulse](http://www.iredex.com/micropulse)

Figure 6. IQ 577 MicroPulse Treatment of Macular Edema Secondary to BRVO

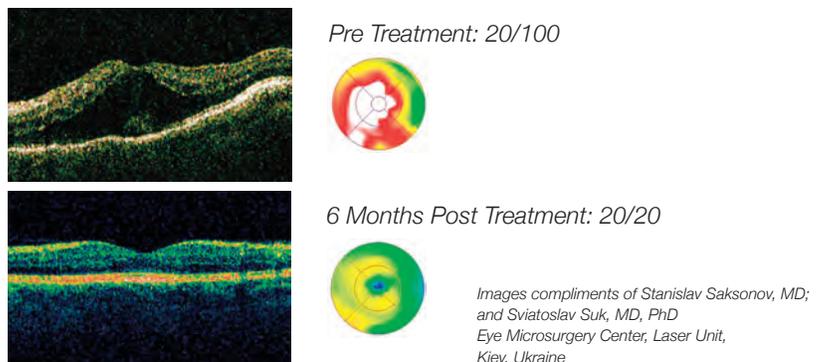

## References

1. L'Esperance FA Jr. Clinical photocoagulation with the organic dye laser. A preliminary communication. Arch Ophthalmol 1985;103(9):1312-6.
2. Mainster MA. Wavelength selection in macular photocoagulation. Tissue optics, thermal effects, and laser systems. Ophthalmology 1986;93(7):952-8
3. Castillejos-Rios D, Devenyl R, Moffat K, Yu E. Dye yellow vs argon green laser in panretinal photocoagulation for proliferative diabetic retinopathy: A comparison of minimum power requirements. Can J Ophthalmol 1992;27(5):243-244
4. Brooks HL, Jr., Eagle RC, Jr., Schroeder RP, Annesley WH, Shields JA, Augsburger JJ. Clinicopathologic study of organic dye. Laser in the human fundus. Ophthalmology 1989;96(6):822-34.
5. Olk RJ, Akduman L. Minimal intensity diode laser (810 nanometer) photocoagulation (MIP) for diffuse diabetic macular edema (ddme). Semin Ophthalmol 2001;16(1):25-30.
6. Moorman CM, Hamilton AM. Clinical applications of the micropulse diode laser. Eye 1999;13 (Pt 2):145-50.
7. Laursen ML, Moeller F, Sander B, Sjoelie AK. Subthreshold micropulse diode laser treatment in diabetic macular oedema. Br J Ophthalmol 2004;88(9):1173-9.
8. Tseng Shih-Yu. Clinical application of micropulse diode laser in the treatment of macular edema. Am J Ophthalmol 2005;139(4):S58.
9. Luttrull JK, Musch DC, Mainster MA. Subthreshold diode micropulse photocoagulation for the treatment of clinically significant diabetic macular oedema. Br J Ophthalmol 2005;89(1):74-80.
10. Luttrull JK, Spink CJ. Serial optical coherence tomography of subthreshold diode laser micropulse photocoagulation for diabetic macular edema. Ophthalmic Surg Lasers Imaging 2006;37(5):370-7.
11. Sivaprasad S, Sandhu R, Tandon A, Sayed-Ahmed K, McHugh DA. Subthreshold micropulse diode laser photocoagulation for clinically significant diabetic macular oedema: A three-year follow up. Clin Experiment Ophthalmol 2007;35(7):640-4.
12. Nakamura Y, Tatsumi T, Arai M, Takatsuna Y, Mitamura Y, Yamamoto S. [subthreshold micropulse diode laser photocoagulation for diabetic macular edema with hard exudates]. Nippon Ganka Gakkai Zasshi 2009;113(8):787-91.
13. Figueira J, Khan J, Nunes S, Sivaprasad S, Rosa A, de Abreu JF, Cunha-Vaz JG, Chong NV. Prospective randomised controlled trial comparing sub-threshold micropulse diode laser photocoagulation and conventional green laser for clinically significant diabetic macular oedema. Br J Ophthalmol 2009;93(10):1341-4.
14. Ohkoshi K, Yamaguchi T. Subthreshold micropulse diode laser photocoagulation for diabetic macular edema in Japanese patients. Am J Ophthalmol 2010;149(1):133-9.
15. Vujosevic S, Bottega E, Casciano M, Pilotto E, Convento E, Midena E. Microperimetry and fundus autofluorescence in diabetic macular edema: Subthreshold micropulse diode laser versus modified early treatment diabetic retinopathy study laser photocoagulation. Retina 2010;30(6):908-916.
16. Kumar V, Ghosh B, Mehta DK, Goel N. Functional outcome of subthreshold versus threshold diode laser photocoagulation in diabetic macular oedema. Eye (Lond) 2010;24(9):1459-65.
17. Lavinsky D, Cardillo JA, Melo LA, Jr., Dare A, Farah ME, Belfort R Jr. Randomized clinical trial evaluating mETDRS versus normal or high-density micropulse photocoagulation for diabetic macular edema. Invest Ophthalmol Vis Sci 52(7):4314-23.
18. Luttrull JK, Sramek C, Palanker D, Spink CJ, Musch DC. Long-term safety, high-resolution imaging, and tissue temperature modeling of sub-visible diode micropulse photocoagulation for retinovascular macular edema. Retina 2011; Publish Ahead of Print:10.1097/IAE.0b013e3182206f6c.
19. Venkatesh P, Ramanjulu R, Azad R, Vohra R, Garg S. Subthreshold micropulse diode laser and double frequency neodymium: YAG laser in treatment of diabetic macular edema: A prospective, randomized study using multifocal electroretinography. Photomed Laser Surg 2011.
20. Luttrull JK, Musch DC, Spink CA. Subthreshold diode micropulse panretinal photocoagulation for proliferative diabetic retinopathy. Eye (Lond) 2008;22(5):607-12.
21. Kumar V, Ghosh B, Raina UK, Goel N. Subthreshold diode micropulse panretinal photocoagulation for proliferative diabetic retinopathy. Eye 2009;23(11):2122-2123.
22. Parodi MB, Spasse S, Iacono P, Di Stefano G, Canziani T, Ravalico G. Subthreshold grid laser treatment of macular edema secondary to branch retinal vein occlusion with micropulse infrared (810 nanometer) diode laser. Ophthalmology 2006;113(12):2237-42.
23. Parodi MB, Iacono P, Ravalico G. Intravitreal triamcinolone acetonide combined with subthreshold grid laser treatment for macular oedema in branch retinal vein occlusion: A pilot study. Br J Ophthalmol 2008;92(8):1046-50.
24. Luttrull JK. Laser for BRVO: History and current practice. Retina Today 2011;May/June:74-76.
25. Ricci F, Missiroli F, Cerulli L. Indocyanine green dye-enhanced micropulsed diode laser: A novel approach to subthreshold RPE treatment in a case of central serous chorioretinopathy. Eur J Ophthalmol 2004;14(1):74-82.
26. Lanzetta P, Furlan F, Morgante L, Veritti D, Bandello F. Nonvisible subthreshold micropulse diode laser (810 nm) treatment of central serous chorioretinopathy. A pilot study. Eur J Ophthalmol 2008;18(6):934-40.
27. Chen SN, Hwang JF, Tseng LF, Lin CJ. Subthreshold diode micropulse photocoagulation for the treatment of chronic central serous chorioretinopathy with juxtafoveal leakage. Ophthalmology 2008;115(12):2229-34.
28. Gupta B, Elagouz M, McHugh D, Chong V, Sivaprasad S. Micropulse diode laser photocoagulation for central serous chorio-retinopathy. Clin Experiment Ophthalmol 2009;37(8):801-5.
29. Ricci F, Missiroli F, Regine F, Grossi M, Dorin G. Indocyanine green enhanced subthreshold diode-laser micropulse photocoagulation treatment of chronic central serous chorioretinopathy. Graefes Arch Clin Exp Ophthalmol 2009;247(5):597-607.
30. Dare AR, Lavinsky D, Magalhaes F, Roisman L, Tognin F, Moreira CE, Cardillo JA. Focal juxtafoveal and grid pattern selective micro-pulsed laser photocoagulation for treatment of chronic central serous chorioretinopathy. Invest Ophthalmol Vis Sci 2009;50(5):ARVO E-Abstract 214.
31. Maia AM, Penha FM, Regatieri CVS, Cardillo JA, Farah ME. Micropulse 577 nm - yellow laser photocoagulation for central serous chorio-retinopathy. Invest Ophthalmol Vis Sci 2010;51(5):4273-.
32. Dare AJ, Cardillo JA, Lavinsky D, Belfort Rubens Jr., Moreira CE. 577 nm yellow selective subthreshold laser photocoagulation for the treatment of central serous chorioretinopathy with foveal leakage. Invest Ophthalmol Vis Sci 2011;52(6):6622-.
33. Peroni R, Cardillo JA, Dare AJ, Aguirre JG, Lavinsky D, Farah ME, Belfort R Jr. A combined low energy, short pulsed 577 nm mild macular grid photocoagulation with 577 nm-micropulsed central laser stimulation for diabetic macular edema with foveal leakage (the sandwich grid). Invest Ophthalmol Vis Sci 2011;52(6):590-.
34. Aguirre JGM, Sr., Cardillo JA, Dare AJ, Peroni R, Lavinsky D, Farah ME, Belfort Rubens Jr. 577 nm short pulsed and low energy selective macular grid laser photocoagulation for diffuse diabetic macular edema. Invest Ophthalmol Vis Sci 2011;52(6):592-.
35. Pasechnikova N, Suk S. 577 nm micropulse laser treatment of macular edema secondary to branch retinal vein occlusion. Retina Today 2011;Supplement (April):11-13.
36. Mainster MA, Whitacre MM. Dye yellow photocoagulation of retinal arterial macroaneurysms. Am J Ophthalmol 1988;105(1):97-8.

 Emergo Europe  
Molenstraat 15, 2513 BH The Hague, The Netherlands, Tel.: (31) (0) 70 345-8570, Fax: (31) (0) 70 346-7299

U.S. Patent No. 7,771,417

